


PRODUCT INFORMATION  
**IL-17D (Interleukin-17D), Human**

<b>Catalog number</b>	C01021-5UG / C01021-20UG / C01021-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	The protein encoded by this gene is a cytokine that shares the sequence similarity with IL17. The treatment of endothelial cells with this cytokine has been shown to stimulate the production of other cytokines including IL6, IL8 and CSF2/ GM-CSF. The increased expression of IL8 induced by this cytokine was found to be NF-kappa B-dependent.
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	APRAGRPARPRGCADRPEELLEQLYGRLAAGVLSAFHHTLQLGPREQARNA SCPAGGRPADDRFRPPTNLRVSPWAYRISYDPARYPRYLPEAYCLCRGCLTG LFGEDVRFRSAPVYMPTVLRRTACAGGRSVYTEAYVTIPVGCTCVPEPEK DADSINSSIDKQGAKLLLGPNDAPAGP with polyhistidine tag at the N-terminus
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Purity</b>	>98% as determined by SDS-PAGE. Ni-NTA chromatography
<b>Formulation</b>	The protein was lyophilized from a solution containing 1X PBS, pH 8.0.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Storage</b>	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
<b>Note</b>	Please use within two weeks after protein reconstitution.

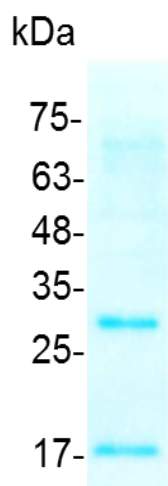

SDS-PAGE analysis of recombinant human IL-17D

For Research Use Only.