


PRODUCT INFORMATION  
**IL-29 (Interleukin-29), Human**

<b>Catalog number</b>	C01035-5UG / C01035-20UG / C01035-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	Interleukin-29 (IL-29) is a member of the helical cytokine family and is a type III interferon. It is also known as IFNλ1 and is highly similar in amino acid sequence to the IL-28, the other type III interferon. IL-29 plays an important role in host defenses against microbes and its gene is highly upregulated in cells infected with viruses. IL29 is not present in the mouse genome.
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	MPTSKPTTTGKGCHIGRFKSLSPQELASFKKARDALEESLKLKNWSCSSPVFP GNWDLRLLQVRERPVALEAELALTLKVLEAAAGPALEDVLDQPLHTLHHILSQL QACIQPQPTAGPRPRGRLHHWLHRLQEAPKKESAGCLEASVTFNLFRLLTRDL KYVADGNLCLRTSTHPEST with polyhistidine tag at the C-terminus
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Activity</b>	Measure by its ability to induce IL-8 secretion in HuH7 cells. The ED <sub>50</sub> for this effect is <6 ng/mL.
<b>Purity</b>	>98% as determined by SDS-PAGE. Ni-NTA chromatography
<b>Formulation</b>	The protein was lyophilized from a solution containing 1X PBS, pH 8.0.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Storage</b>	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
<b>Note</b>	Please use within one month after protein reconstitution.

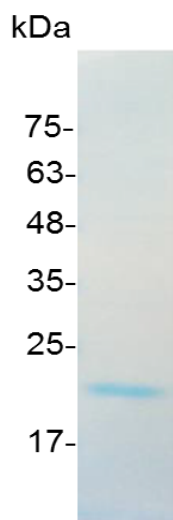

SDS-PAGE analysis of recombinant human IL-29

For Research Use Only.