


PRODUCT INFORMATION

BMP-9 (Bone morphogenetic protein-9), Human

Catalog number	C01070-5UG / C01070-20UG / C01070-100UG
Package	5 µg / 20 µg / 100 µg
Description	BMP-9 is a member of the BMP subgroup of the TGF-beta superfamily proteins that signal through heterodimeric complexes composed of type I and type II BMP receptors. BMP-9 regulates the development and function of a variety of embryonal and adult tissues
Source	<i>Escherichia coli</i>
Sequence	SAGAGSHCQKTSLRVNFEDIGWDSWIIAPKEYEAYECKGGCFFPLADDVTPTK HAIVQTLVHLKFPTKVGKACCVPTKLSPIVLYKDDMGVPTLKYHYEGMSVAEC GCR with polyhistidine tag at the N-terminus
Endotoxin level	<0.01 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to induce alkaline phosphatase production by ATDC5 cells. The ED ₅₀ for this effect is <0.4 ng/mL.
Purity	>98% as determined by SDS-PAGE. Ni-NTA chromatography
Formulation	The protein was lyophilized from a solution containing 20 mM sodium citrate, 0.2 M NaCl, pH 3.5.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

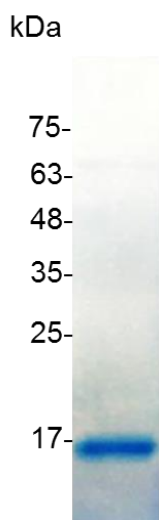

SDS-PAGE analysis of recombinant human BMP-9

For Research Use Only.