


## PRODUCT INFORMATION

**IL-12 p35 (Interleukin-12 p35), Mouse**

<b>Catalog number</b>	C02013-5UG / C02013-20UG / C02013-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	Active IL-12 is a p70 disulfide-linked dimer composed of p35 and p40 subunits. The protein is a pleiotropic cytokine produced primarily by antigen presenting cells and has multiple effects on T lymphocytes and natural killer cells in terms of stimulating cytotoxicity, proliferation, production of other cytokines and Th1 subset differentiation.
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	MRVIPVSGPARCLSQSRNLLKTTDDMVKTAREKCLKHYSCTAEDIDHEDITRDQT STLKTCLPLELHKNESCLATRETSSTTRGSLPPQKTSMMTLCLGSIYEDLKM YQTEFQAINAALQNHNHQQIILDKGMLVAIDELMQSLNHNGETLRQKPPVGEAD PYRVKMKLCILLHAFSTRVVTINRVMGYLSSA with polyhistidine tag at the C-terminus
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Activity</b>	Measure by its ability to induce proliferation in T-cell enriched PBMC. The ED <sub>50</sub> for this effect is <0.2 ng/mL.
<b>Purity</b>	>98% as determined by SDS-PAGE. Ni-NTA chromatography
<b>Formulation</b>	The protein was lyophilized from a solution containing 1X PBS, pH 7.4.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Storage</b>	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
<b>Note</b>	Please use within one month after protein reconstitution.

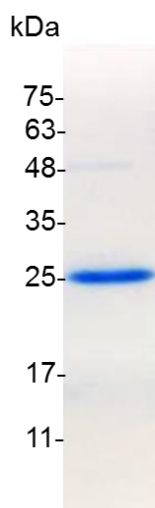

SDS-PAGE analysis of recombinant mouse IL-12 p35

For Research Use Only.