


PRODUCT INFORMATION
IL-17B (Interleukin-17B), Mouse

Catalog number	C02019-5UG / C02019-20UG / C02019-100UG
Package	5 µg / 20 µg / 100 µg
Description	The interleukin (IL)-17 superfamily, a relatively new family of cytokines, consists of six ligands (from IL-17A to IL-17F), which bind to five receptor subtypes (from IL-17RA to IL-17RE) and induce downstream signaling. L17B Diseases associated with IL17B include spondyloarthritis and neuronitis, and among its related super-pathways are Mucin expression in CF via IL-6, IL-17 signaling pathways, and STAT3 Pathway.
Source	<i>Escherichia coli</i>
Sequence	MHPRNTKGKRKGQGRPSPLAPGPHQVPLDLVSRVKPYARMEEYERNLGEMV AQLRNSSEPAKKKCEVNLQLWLSNKRSLSPWGYSHDPSRIPADLPEARCLCL GCVNPFTMQEDRSMVSPVFSQVPVRRRLCPQPPRPGPCRQRVVMETIAVGC TCIF with polyhistidine tag at the C-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to induce IL-8 secretion in HepG2 cells. The ED ₅₀ for this effect is <1.5 ng/mL.
Purity	>98% as determined by SDS-PAGE. Ni-NTA chromatography
Formulation	The protein was lyophilized from a solution containing 1X PBS, pH 7.4.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

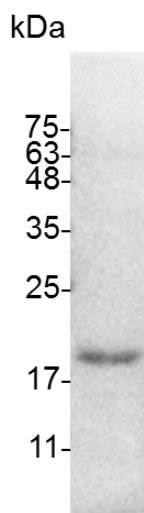

SDS-PAGE analysis of recombinant mouseIL-17B

For Research Use Only.