


## PRODUCT INFORMATION

**beta-NGF (Nerve growth factor-beta), Mouse**

<b>Catalog number</b>	C02092-5UG / C02092-20UG / C02092-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	Nerve growth factor (NGF) is a neurotrophic factor and neuropeptide primarily involved in the regulation of growth, maintenance, proliferation, and survival of certain target neurons. NGF-β acts through its receptor β-NGFR and is involved in the development and maintenance of the sensory and sympathetic nervous systems. NGF-β also is also involved in the growth, differentiation, and survival of B lymphocytes. Human, mouse and rat proteins show cross-reactivity.
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	MSSTHPVFHMGFEFSVCDSVSVVVGDKTTATDIKGKEVTVLAEVNINNSVFRQY FFETKCRASNPVESGCRGIDSKHWNSYCTTTHTFVKALTTDEKQAAWRFRIDT ACVCVLSRKATRRG with polyhistidine tag at the C-terminus
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Activity</b>	Measure by its ability to induce TF-1 cells proliferation. The ED <sub>50</sub> for this effect is <1 ng/mL. The specific activity of recombinant mouse beta-NGF is > 1 x 10 <sup>6</sup> IU/mg.
<b>Purity</b>	>98% as determined by SDS-PAGE. Ni-NTA chromatography
<b>Formulation</b>	The protein was lyophilized from a solution containing 20 mM sodium citrate, 0.2 M NaCl, pH 4.5.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Storage</b>	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
<b>Note</b>	Please use within one month after protein reconstitution.

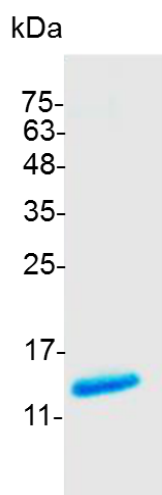

SDS-PAGE analysis of recombinant mouse beta-NGF

For Research Use Only.