


PRODUCT INFORMATION
IL-2 (Interleukin-2), Swine

Catalog number	C03003-5UG / C03003-20UG / C03003-100UG
Package	5 µg / 20 µg / 100 µg
Description	Interleukin-2 (IL-2) is an interleukin, a type of cytokine signaling molecule in the immune system. It is a 15,5 - 16 kDa protein that regulates the activities of white blood cells (leukocytes, often lymphocytes) that are responsible for immunity. IL-2 is part of the body's natural response to microbial infection, and in discriminating between foreign ("non-self") and "self". IL-2 mediates its effects by binding to IL-2 receptors, which are expressed by lymphocytes.
Source	<i>Escherichia coli</i>
Sequence	MAPTSSSTKNTKKQLEPLLLDLQLLLKEVKNYENADLSRMLTFKIFYMPKQATEL KHLQCLVEELKALEGVLNLGQSKNSDSANIKESMNNINVTVLELKGSETSFKCE YDDETVAVEFLNKWITFCQSIYSTLT with polyhistidine tag at the C-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to induce proliferation in CTLL2 cells. The ED ₅₀ for this effect is <0.5 ng/mL.
Purity	>98% as determined by SDS-PAGE. Ni-NTA chromatography
Formulation	The protein was lyophilized from a solution containing 1X PBS, pH 7.4.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

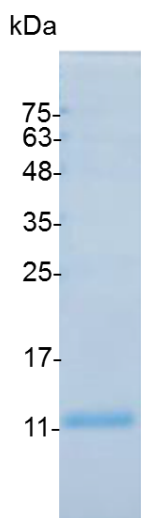

SDS-PAGE analysis of recombinant swine IL-2

For Research Use Only.