


## PRODUCT INFORMATION

**BMP-4 (Bone morphogenetic protein-4), Swine**

<b>Catalog number</b>	C03011-5UG / C03011-20UG / C03011-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	Bone morphogenetic protein 4 is a protein that in humans is encoded by BMP4 gene. BMP4 is found on chromosome 14q22-q23. BMP4 is a member of the bone morphogenetic protein family which is part of the transforming growth factor-beta superfamily. The superfamily includes large families of growth and differentiation factors. BMP4 is highly conserved evolutionarily. BMP4 is found in early embryonic development in the ventral marginal zone and in the eye, heart blood and otic vesicle.
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	MREKQPNYGLAIEVTHLHQTRTHQGQHVRIISRSLPQGGSDWAQLRPLLVTFGH DGRGHALTRRRRAKRSPKHHPQRARKKKNKNCRRHSLYVDFSDVGWNDWIVA PPGYQAFYCHGDCPFPLADHLNSTNHAIVQTLVNSVNSSIPKACCVPTLSAIS MLYLDEYDKVVLKKNYQEMVVEGCGC with polyhistidine tag at the C-terminus
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Purity</b>	>98% as determined by SDS-PAGE. Ni-NTA chromatography
<b>Formulation</b>	The protein was lyophilized from a solution containing 20 mM sodium citrate, 0.2 M NaCl pH 3.5.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Storage</b>	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
<b>Note</b>	Please use within one month after protein reconstitution.

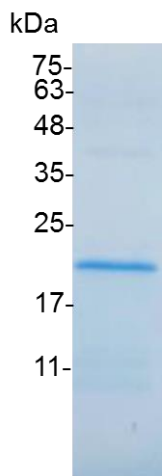

SDS-PAGE analysis of recombinant swine BMP-4

For Research Use Only.