


PRODUCT INFORMATION

IFN alpha 1a (Interferon alpha 1a), Human

Catalog number	C01078-5UG / C01078-20UG / C01078-100UG
Package	5 µg / 20 µg / 100 µg
Description	Interferon alpha-1 or interferon alpha-D, encoded by the human IFNA1 gene, is a member of the type I interferon family. It is expressed in peripheral blood leukocytes and lymphoblastoid cells. There are four highly conserved cysteine residues which form two disulfide bonds. One bond is a necessity for biological activity, which plays a critical role in conducting non-specific resistance against a variety of viral infections, and regulation of cell proliferation and modulates immune responses.
Source	<i>Escherichia coli</i>
Sequence	CDLPETHSLDNRRTLMLLAQMSRISPSSCLMDRHDFGFPQEEFDGNQFQKAPA ISVLHELIQQIFNLFITTKDSSAAWDEDLLDKFCTELYQQQLNDLEACVMQEERVGE TPLMNADSILAVKKYFRRITLYLTEKKYSPCAWEVVRAEIMRSLSLSTNLQERLR RKE with polyhistidine tag at the N-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to inhibit IL-8 secretion in human PBMCs in the presence of LPS. The ED ₅₀ for this effect is <1.12 µg /mL.
Purity	>95% as determined by SDS-PAGE. Ni-NTA chromatography
Formulation	The protein was lyophilized from a solution containing 1X PBS, pH 8.0.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

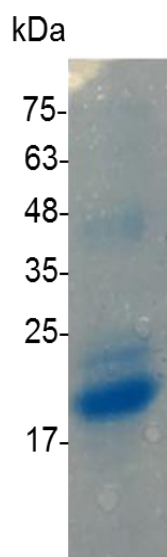

SDS-PAGE analysis of recombinant human IFN alpha 1a

For Research Use Only.