


PRODUCT INFORMATION

CXCL5 (C-X-C motif chemokine 5), Human

Catalog number	C01131-5UG / C01131-20UG / C01131-100UG
Package	5 µg / 20 µg / 100 µg
Description	The protein encoded by this gene, Chemokine (C-X-C motif) ligand 5 (CXCL5), is a small cytokine belonging to the CXC chemokine family that is also known as epithelial-derived neutrophil-activating peptide 78 (ENA-78). This chemokine is produced concomitantly with interleukin-8 (IL8) in response to stimulation with either interleukin-1 (IL1) or tumor necrosis factor-alpha (TNFA). It is observed that, CXCL5 also expresses in eosinophils, and can interact with the type II interferon IFN- γ , thereby cause an inhibition. This chemokine stimulates the chemotaxis of neutrophils possesses angiogenic properties, and elicits these effects by interacting with the cell surface chemokine receptor CXCR2.
Source	<i>Escherichia coli</i>
Sequence	RELRCVCLQTTQG VHPKMISNLQVFAIGPQCSKVEVVASLKNNGKEICLDPEAPF LKKVIQKILDGGNKEN with polyhistidine tag at the N-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to chemoattract BaF3 cells transfected with human CXCR2. The ED ₅₀ for this effect is <10 ng/mL.
Purity	>98% as determined by SDS-PAGE. Ni-NTA chromatography
Formulation	The protein was lyophilized from a solution containing 1X PBS, pH 7.4.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

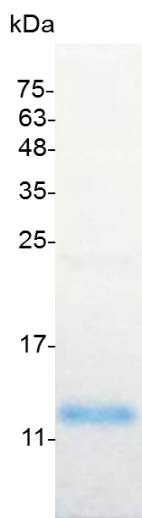

SDS-PAGE analysis of recombinant human CXCL5

For Research Use Only.