


TGF alpha (Transforming growth factor alpha), Human

v. 231001

Catalog number	C01087-5UG / C01087-20UG / C01087-100UG
Package	5 µg / 20 µg / 100 µg
Description	<p>Transforming growth factor alpha (TGF-α) is a protein that in humans is encoded by the TGFA gene. TGF-alpha is a member of the EGF family of cytokines that are synthesized as transmembrane precursors and are characterized by the presence of one or several EGF structural units in their extracellular domain. Expression of TGF-alpha is widespread in tumors and transformed cells. TGF-alpha is also expressed in normal tissues during embryogenesis and in adult tissues, including pituitary, brain, keratinocytes, and macrophages.</p>
Source	<i>Escherichia coli</i>
Sequence	MVVSFNDPCDSHTQFCFHGTCRFLVQEDKPACVCHSGYVGARCEHADLLA with polyhistidine tag at the C-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to induce 3T3 cells proliferation. The ED ₅₀ for this effect is <0.2 ng/mL. The specific activity of recombinant human TGF alpha is > 5 x 10 ⁶ IU/mg.
Purity	>98% as determined by SDS-PAGE.
Form	Lyophilized
Storage Buffer	Lyophilized from a 0.2 µm filtered solution of PBS, pH 8.0.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 200 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Stability & Storage	<p>This product is stable after storage at:</p> <ul style="list-style-type: none"> -20°C for 12 months in lyophilized state from date of receipt. -20°C or -80°C for 1 month under sterile conditions after reconstitution. <p>Avoid repeated freeze/thaw cycles.</p>

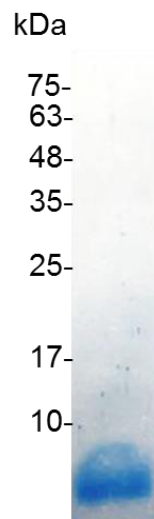

SDS-PAGE analysis of recombinant human TGF alpha

For research use only.