


**IL-17F (Interleukin-17F), Mouse**

v. 231101

<b>Catalog number</b>	C02021-5UG / C02021-20UG / C02021-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	The protein encoded by this gene is a cytokine that shares sequence similarity with IL17. This cytokine is expressed by activated T cells, and has been shown to stimulate the production of several other cytokines, including IL6, IL8, and CSF2/GM-CSF. This cytokine is also found to inhibit the angiogenesis of endothelial cells and induce endothelial cells to produce IL2, TGFB1/TGFB, and monocyte chemoattractant protein-1.
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	MRKNPKAGVPALQKAGNCPPLEDNTVRVDIRIFNQNGGISVPREFQNRSSSPW DYNITRDPHRFPSEIAEAQCRHSGCINAQQQEDSTMNSVAIQQEILVLRREPQG CSNSFRLEKMLLKVGCTCVKPIVHQAA with polyhistidine tag at the C-terminus
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Activity</b>	Measure by its ability to induce IL-6 secretion in 3T3 cells. The ED <sub>50</sub> for this effect is <100 ng/mL.
<b>Purity</b>	>98% as determined by SDS-PAGE.
<b>Form</b>	Lyophilized
<b>Storage Buffer</b>	Lyophilized from a 0.2 µm filtered solution containing 20 mM sodium citrate, pH 4.5.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 200 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Stability &amp; Storage</b>	This product is stable after storage at: <ul style="list-style-type: none"> <li>-20°C for 12 months in lyophilized state from date of receipt.</li> <li>-20°C or -80°C for 1 month under sterile conditions after reconstitution.</li> </ul> Avoid repeated freeze/thaw cycles.

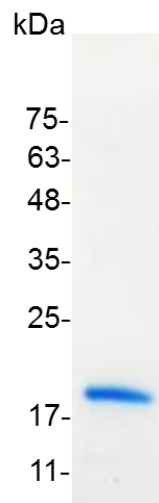

SDS-PAGE analysis of recombinant mouse IL-17F

*For Research Use Only.*