


**VEGF121 (Vascular endothelial growth factor 121), Mouse**

v. 231101

<b>Catalog number</b>	C02071-5UG / C2071-20UG / C02071-100UG
<b>Package</b>	5 µg / 20 µg / 100 µg
<b>Description</b>	<p>Vascular endothelial growth factor (VEGF), originally known as vascular permeability factor (VPF), is a signal protein produced by cells that stimulates the formation of blood vessels. VEGF is required during embryogenesis to regulate the proliferation, migration, and survival of endothelial cells. In adults, VEGF functions mainly in wound healing and the female reproductive cycle. Pathologically, it is involved in tumor angiogenesis and vascular leakage. Circulating VEGF levels correlate with disease activity in autoimmune diseases such as rheumatoid arthritis, multiple sclerosis and systemic lupus erythematosus. VEGF is induced by hypoxia and cytokines such as IL-1, IL-6, IL-8, oncostatin M and TNF-alpha.</p>
<b>Source</b>	<i>Escherichia coli</i>
<b>Sequence</b>	<p>MAPTTEGEQKSHEVIKFMVDVYQRSYCRPIETLVDIFQEYPDEIEYIFKPSCVPLM          RCAGCCNDEALECVPTSESNITMQIMRIKPHQSQHIGEMSFLQHSRCECRPKK          DRTKPEKCDKPRR with polyhistidine tag at the C-terminus</p>
<b>Endotoxin level</b>	<0.1 EU per 1 µg of the protein by the LAL method.
<b>Activity</b>	Measure by its ability to induce proliferation in HUVEC cells. The ED <sub>50</sub> for this effect is <3 ng/mL.
<b>Purity</b>	>95% as determined by SDS-PAGE.
<b>Form</b>	Lyophilized
<b>Storage Buffer</b>	Lyophilized from a 0.2 µm filtered solution of PBS, pH 8.0.
<b>Reconstitution</b>	It is recommended to reconstitute the lyophilized protein in sterile H <sub>2</sub> O to a concentration not less than 200 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
<b>Stability &amp; Storage</b>	<p>This product is stable after storage at:</p> <ul style="list-style-type: none"> <li>-20°C for 12 months in lyophilized state from date of receipt.</li> <li>-20°C or -80°C for 1 month under sterile conditions after reconstitution.</li> </ul> <p>Avoid repeated freeze/thaw cycles.</p>

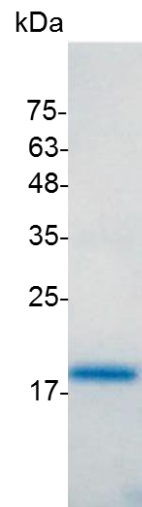

SDS-PAGE analysis of recombinant mouse VEGF121

---

*For Research Use Only.*