


CXCL5 (C-X-C motif chemokine 5), Mouse

v. 231101

Catalog number	C02076-5UG / C02076-20UG / C02076-100UG
Package	5 µg / 20 µg / 100 µg
Description	<p>The protein encoded by this gene, Chemokine (C-X-C motif) ligand 5 (CXCL5), is a small cytokine belonging to the CXC chemokine family that is also known as epithelial-derived neutrophil-activating peptide 78 (ENA-78). This chemokine is produced concomitantly with interleukin-8 (IL8) in response to stimulation with either interleukin-1 (IL1) or tumor necrosis factor-alpha (TNFA). It is observed that, CXCL5 also expresses in eosinophils, and can interact with the type II interferon IFN-γ, thereby cause an inhibition. This chemokine stimulates the chemotaxis of neutrophils possesses angiogenic properties, and elicits these effects by interacting with the cell surface chemokine receptor CXCR2.</p>
Source	<i>Escherichia coli</i>
Sequence	VIAATELRCVCLTVTPKINPKLIANLEVIPAGPQCPTVEVIAKLNQKEVCLDPEA PVIKKIIQKILGSDKKKA with polyhistidine tag at the N-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to chemoattract BaF3 cells transfected with human CXCR2. The ED ₅₀ for this effect is <100 ng/mL.
Purity	>98% as determined by SDS-PAGE.
Form	Lyophilized
Storage Buffer	Lyophilized from a 0.2 µm filtered solution of PBS, pH 7.4.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 200 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Stability & Storage	<p>This product is stable after storage at:</p> <ul style="list-style-type: none"> -20°C for 12 months in lyophilized state from date of receipt. -20°C or -80°C for 1 month under sterile conditions after reconstitution. <p>Avoid repeated freeze/thaw cycles.</p>

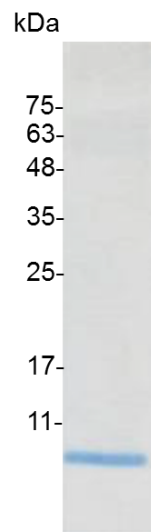

SDS-PAGE analysis of recombinant mouse CXCL5

For Research Use Only.